
How to Capture Multiple Images?

Droplet Lab Projects | Surface Science Instruments & Solutions · https://dropletlab.com/knowledge-doc/how-to-capture-multiple-images/


1. Once you have created a project, choose that project by tapping the project name from the Projects screen.




2. As the project opens, tap the camera icon located at the bottom-left corner of the screen. Doing so opens the camera for you to capture images.

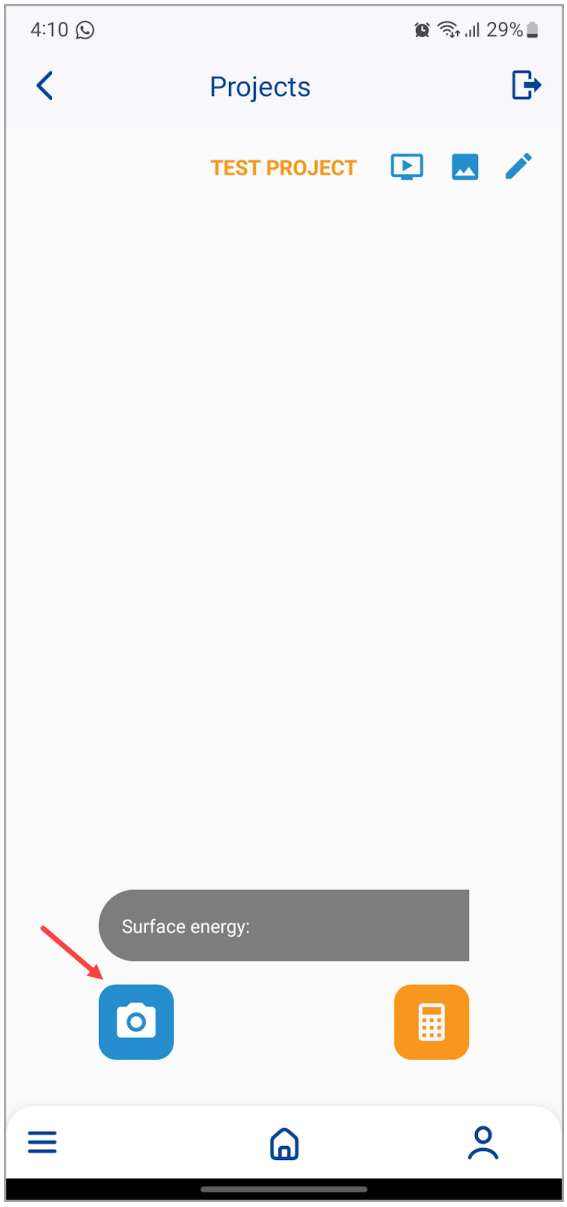


3. As the camera opens, tap the MULTI icon located at the top-right corner of the screen.

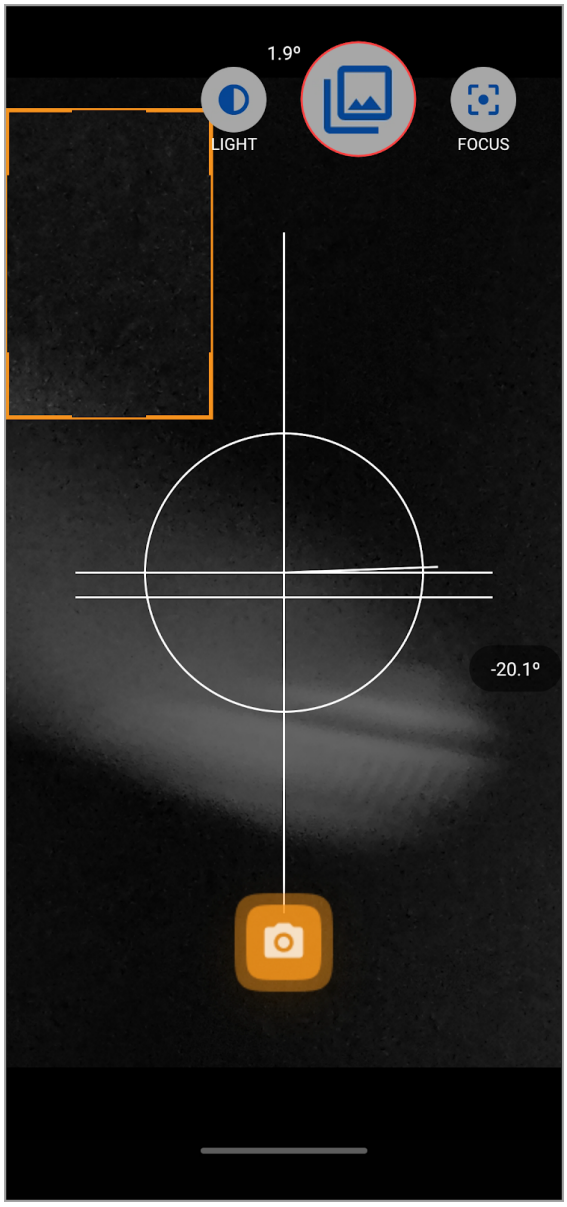


4. As a result, the Multiple Image Capture popup opens. Here, use the Images to Capture field to define the number of images you want to capture.

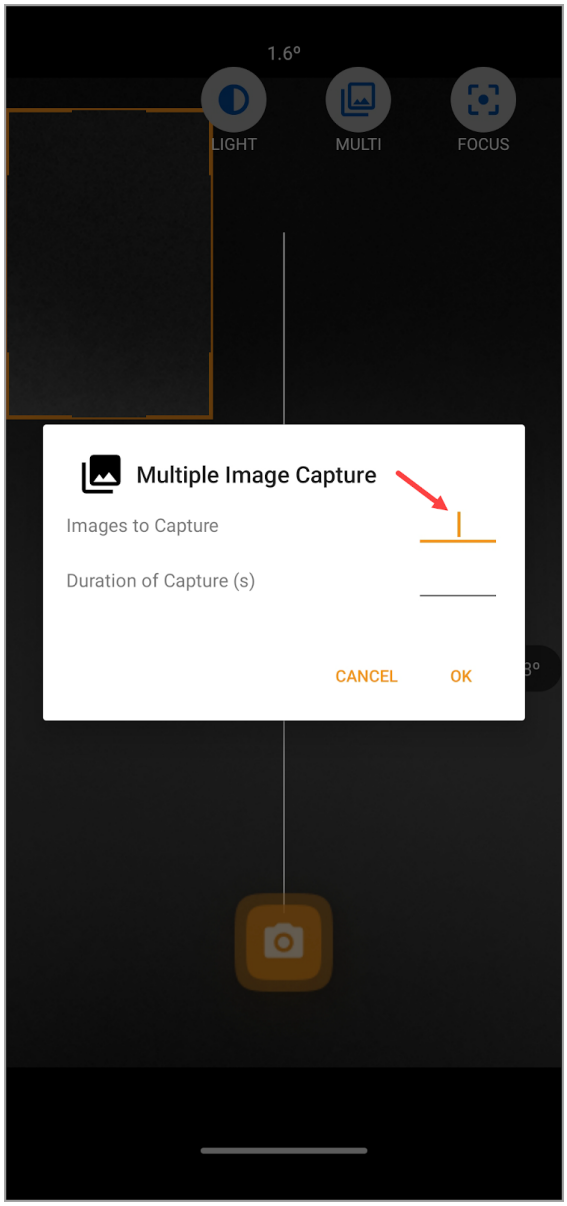

5. Enter the duration for capture(s) using the Duration of Capture(s) field.

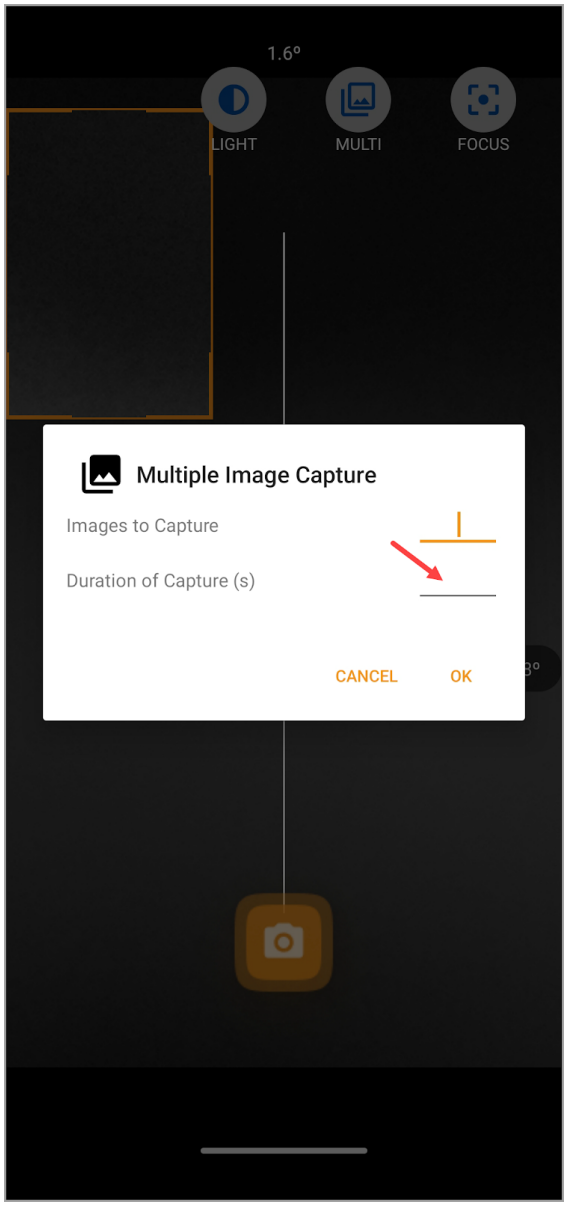

6. Once the required information is filled in, tap the OK button.

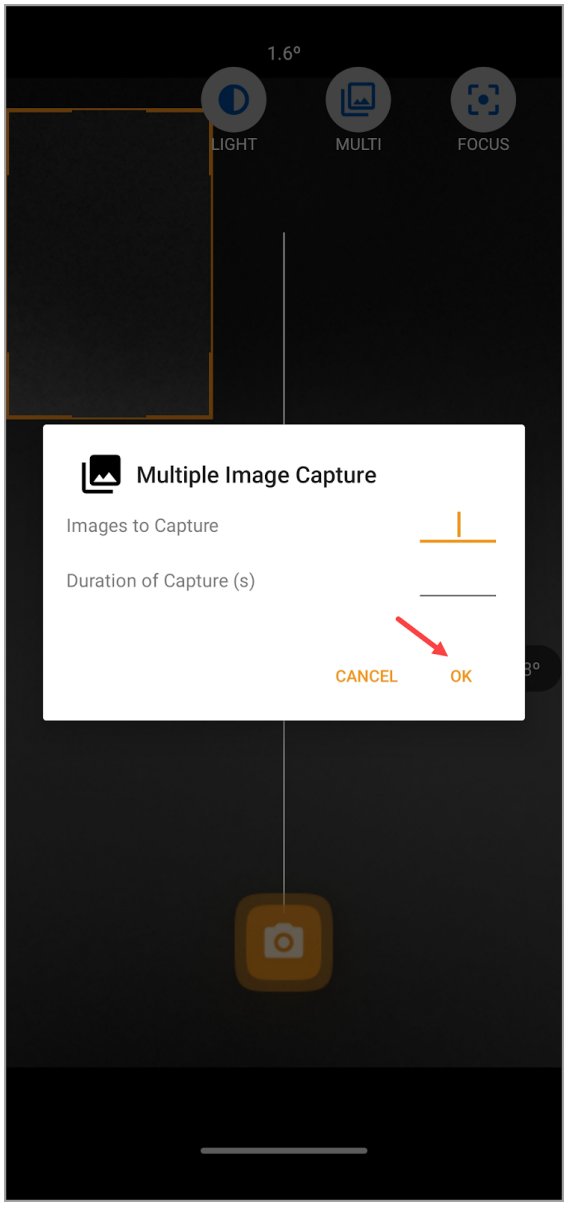
1.6°
LIGHT
MULTI
FOCUS
Multiple Image Capture
Images to Capture
Duration of Capture (s)
CANCEL
OK