
Setting up Application Permissions


App Installation · https://dropletlab.com/knowledge-doc/setting-up-application-permissions/


1. Once the app is downloaded, tap the app icon from the home screen to open the application.

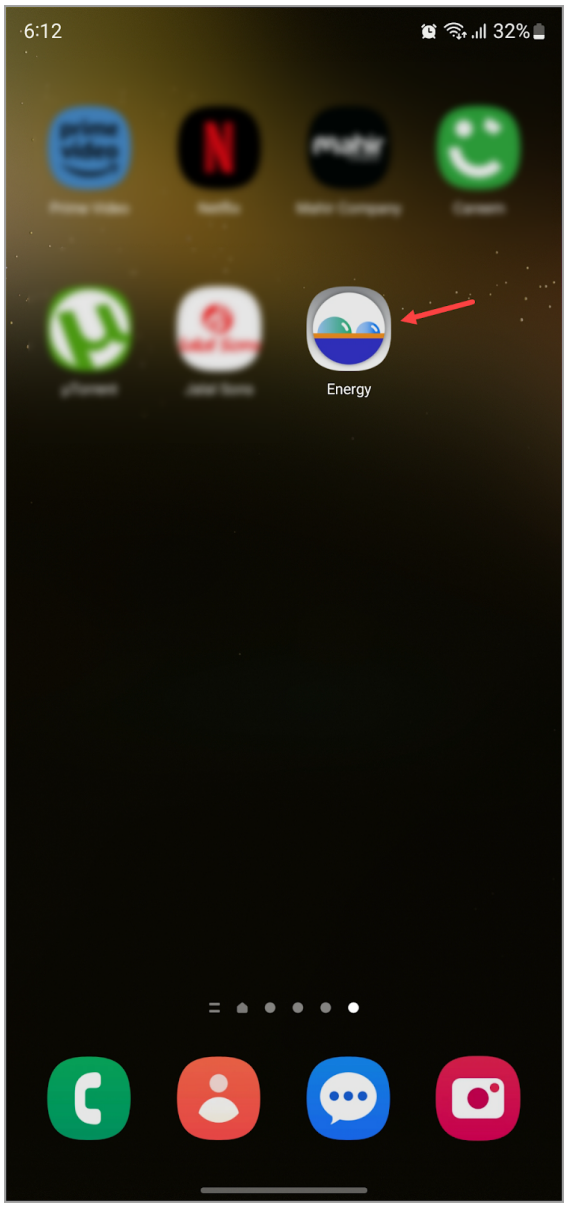


2. As the app opens, it asks for a number of permissions it requires to run. They include the permission to access the device camera so that droplet images can be captured, access system storage so that you can import images, and manage phone calls so that the device’s

unique identifier can be accessed to activate the app. To enable each permission, tap the Allow button.

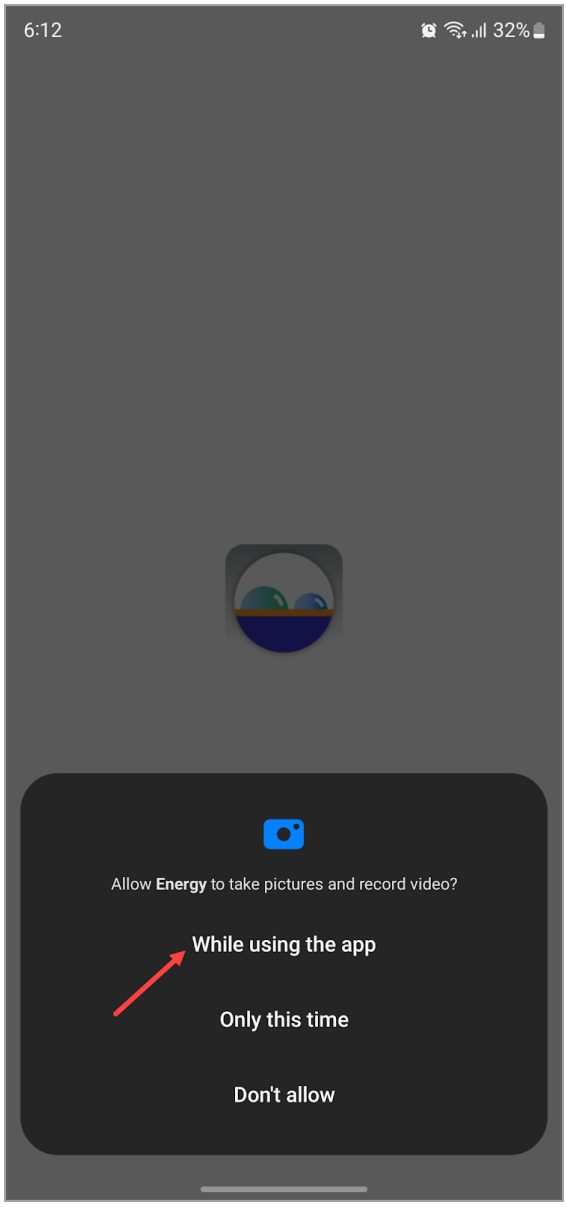

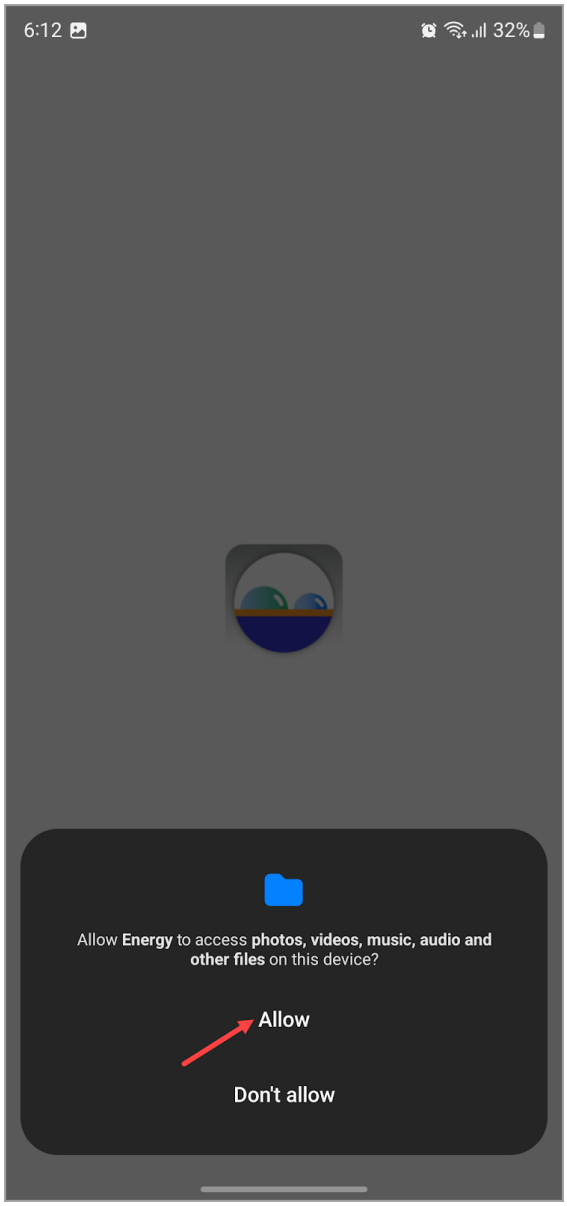
6:12
32%
Allow Energy to access photos, videos, music, audio and other files on this device?
Allow
Don't allow

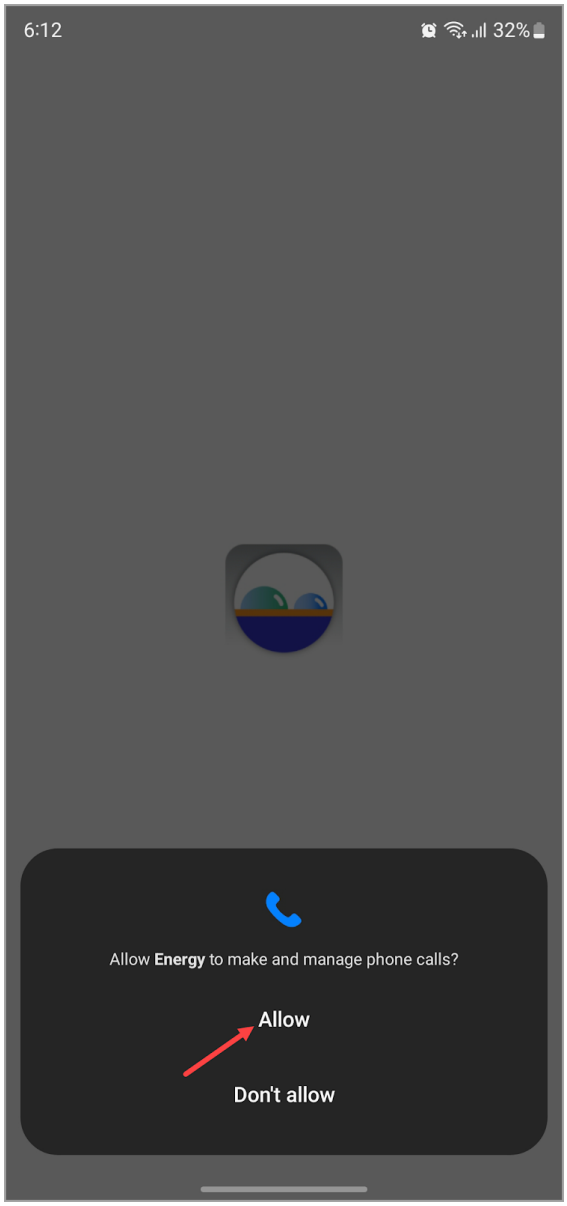

3. As the permissions are allowed, you are taken to the application’s landing screen.